
SINGLE USE HASSON/BLUNT PORT SYSTEM INSTRUCTIONS FOR USE

Including sterile packaged Single-Use Cannula and Hasson Trocar System

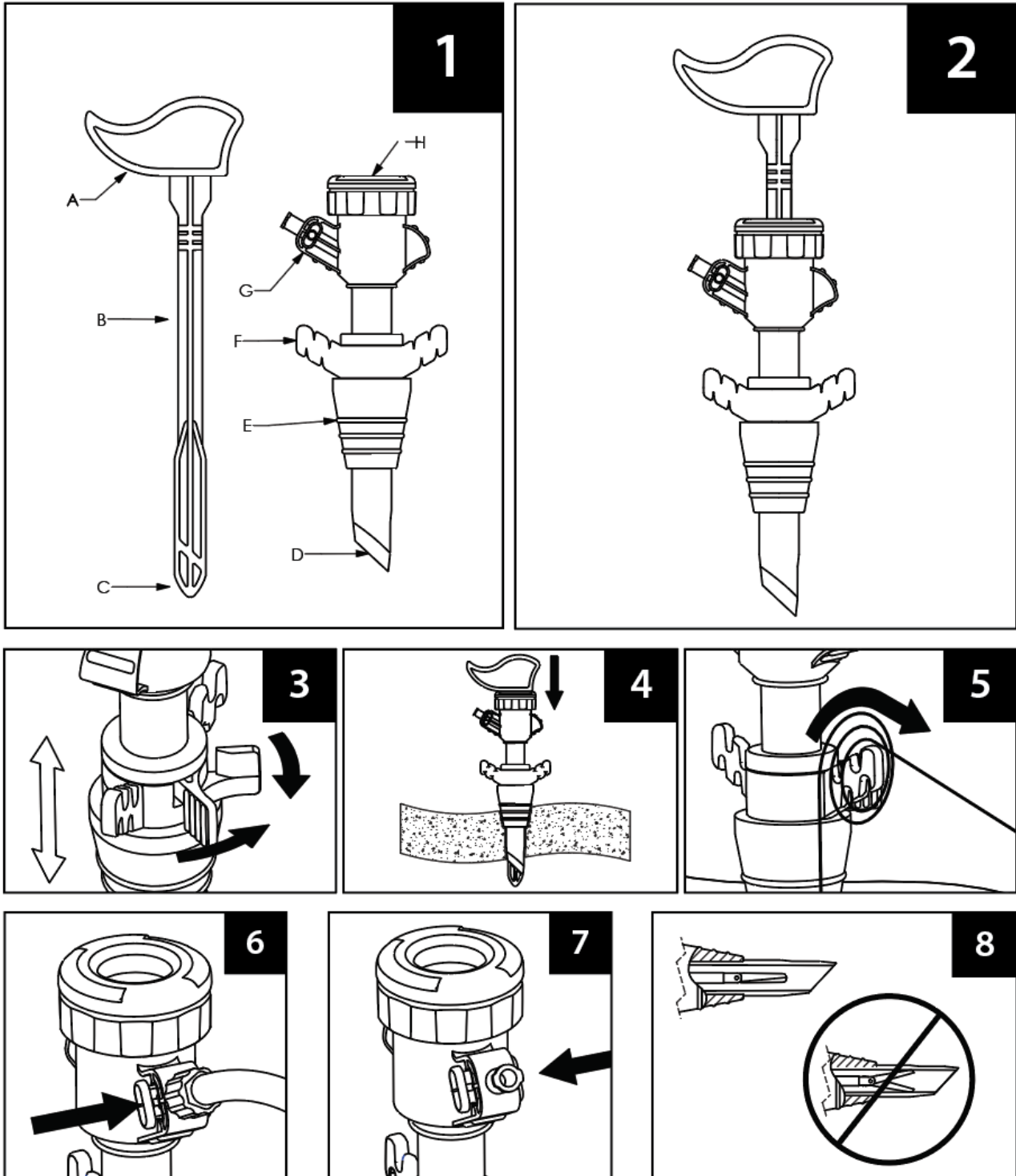BEFORE USING THIS PRODUCT, READ THE FOLLOWING INFORMATION THOROUGHLY.

IMPORTANT!

1. This booklet is designed to assist using this product. It is not a reference to surgical techniques.
2. This device was designed, tested and manufactured for single use patient use only.
3. Reuse of this device may lead to its failure and subsequent patient injury.
4. Reprocessing and/or re-sterilization of this device may create risk of contamination and patient infection.
5. Do not reuse, reprocess or re-sterilize this device

DESCRIPTION/INDICATIONS FOR USE

GENICON trocar/cannula systems are intended for surgical insertion into the peritoneal cavity for the purpose of establishing a port of entry for instrumentation and for visualization during laparoscopic for primary and secondary insertions. These devices have applications in gynecologic, general, thoracic, and urologic endoscopic procedures.

CONTRAINDICATIONS

The use of GENICON laparoscopic and endoscopic products is contraindicated whenever endoscopic surgical techniques are contraindicated for any reason. Contraindications relevant to individual GENICON products are noted in the specific sections.

WARNINGS AND PRECAUTIONS

1. Endoscopic surgery should be performed only by physicians who are thoroughly trained in endoscopic techniques and failure modes, precautions, and corrective actions in the event of failure.
2. GENICON single-use products are intended for single patient use only – DO NOT RESTERILIZE.
3. Do not use any pre-sterilized products if the packaging sterile barrier is damaged
4. Consult medical literature or country specific regulations for specific techniques, complications, and hazards prior to procedure.
5. Care must be taken when using laparoscopic instrumentation to avoid damage to major vessels and other anatomic structures.
6. Establish and maintain adequate pneumoperitoneum to reduce the risk of injury to internal structures.
7. Properly position the patient and note anatomical landmarks to introduced devices without unintended damage.
8. Do not use excessive force or in a manner not consistent with normal instrumentation use.
9. A thorough understanding of the principles and techniques involved in laparoscopic electrosurgical procedures is essential to avoid shock and burn hazards to both patient and medical personnel, and damage to the device or other surgical instruments.
10. Carry out electrosurgical procedures using only non-combustible gases (CO₂).
11. Verify that the devices are compatible with other products that will be used in surgery prior to the procedure.
12. Store all products safely in a climate controlled environment and handle with care.
13. All GENICON devices are latex free and sterilized using gamma sterilization, except as noted.
14. Dispose of all used or damaged products using standard hospital practices for biohazard control.
15. Care must be taken when using any trocar/cannula to avoid damage to major vessels and other anatomic structures.
16. The trocar/cannula should be placed under direct endoscopic vision.
17. All instruments should be aligned with the central axis of the cannula when inserting/removing through the seal to minimize the potential of inadvertent damage to the cannula seal. Extra care should be used when inserting and removing angular and/or asymmetrical instruments, such as “J” hooks, Covidien Endo GIA, and clip appliers.

To minimize the risk of such injury, be sure to:

1. Establish and maintain adequate pneumoperitoneum
2. Properly position the patient to help displace organs out of the area of penetration
3. Note important anatomical landmarks
4. Direct the trocar tip away from major vessels and organs
5. Do not use excessive force during insertion
6. Do not use other manufacturer's trocars or cannulae in conjunction with GENICON trocars or cannulae.

POTENTIAL COMPLICATIONS FOR ALL GENICON® PRODUCTS

Potential complications associated with the use of laparoscopic devices include but are not limited to organ or vessel damage/perforation, hemorrhage, hematoma, trocar site hernia, and septicemia/infection.

INSTRUCTIONS FOR USE

1. Inspect all components carefully for damage during shipment, handling or after maintenance. **Refer to picture - -**
 - A. Trocar Handle
 - B. Trocar Shaft
 - C. Blunt Tip
 - D. Cannula Tip
 - E. Suture Cone
 - F. Suture Wings
 - G. CO₂ Luer Port/Valve
 - H. Seal/Cap Assembly
2. Inspect sterile packaging, and verify that sterile barriers are intact. Do not use if package is compromised.
3. Verify that the sizes of all surgical components selected are compatible.
4. Open package using sterile technique and place trocar and cannula (if applicable) on sterile field.
5. Using proper laparoscopy technique, select entry site to ensure the absence of vessels.
6. Assemble by gently pressing the blunt tip obturator through the center of the primary valve assembly in the cannula; applying sterile saline to cannula valve will aid in assembly of trocar into the cannula. **Refer to picture - -**
7. Adjust the cannula anchoring device by squeezing the suture anchor wing tabs together and sliding the cone as needed. **Refer to picture - -**
8. Make an incision large enough to accommodate the size of the trocar sleeve at the site of placement. Bluntly dissect through the fascia and peritoneum using Hasson/open laparoscopy technique. **Refer to picture - -**
9. Using two sutures of adequate tensile strength, pass one through each fascial edge, and tag them.
10. Immediately after insertion of the device, inspect the puncture site for excessive bleeding and assess for organ/vessel damage. Institute appropriate treatment if injury occurs.
11. Remove the obturator. Tighten the anchoring device against the body wall if necessary to secure the system during use.
12. Wrap fascial sutures around the notches on the suture wings at least three times to secure the assembly. **Refer to picture - -**

13. The insufflation (CO₂) valve is packaged in the CLOSED position (red slide valve fully depressed). To insufflate the abdomen, attach the CO₂ supply to the luer connection, and press the green colored slide valve fully. **Refer to picture -** To stop insufflation and close the valve, press the red colored slide valve fully. **Refer to picture -**
14. Instruments should not be opened while inside the cannula. **Refer to picture -**
15. After surgery, remove the suture from the anchoring device. Squeeze the suture anchor wing tabs together, and turn the suture anchor until it loosens from the fascia. **Refer to picture -** Remove the cannula. Use suture to close the trocar site.

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This product may be covered by one or more U.S. patents and their foreign counter parts. Other Patents may be pending. GENICON and the GENICON logo design are U.S. registered trademarks of GENICON Rights Reserved.

Disclaimer of Warranties on Re-Processing or Re-Sterilizing of Vendor's Disposable (Single-Patient) Products: Any reprocessing or re-sterilizing of disposable single patient devices voids Vendor's warranty. Vendor assumes no liability for reprocessed or re-sterilized products. Vendor warrants that reasonable care has been used in the manufacture of these disposable single patient devices and that they are free from defects in workmanship or materials at the time of shipment from Vendor. Vendor's sole obligation shall therefore be to repair or replace any device which it determines was defective at the time of shipment. Because no product is completely effective under all circumstances, and because the actual use and handling of these devices are beyond Vendor's control, Vendor cannot warrant for a good effect or against a bad effect in the application and use of these disposable single patient devices. The Purchaser therefore assumes all liability resulting from re-sterilization or re-processing of these disposable single patient products. Vendor therefore gives no warranty of merchantability or fitness for a particular purpose. Vendor will not be liable for any direct, indirect, incidental or consequential damages resulting from or related to a disposable product that has been re-sterilized or re-processed. Vendor shall not be liable for incidental or consequential loss, damage or expense resulting from the use or application of these said products. This warranty is in lieu of all other warranties, whether implied, express, oral or written, and no individual or entity has the authority to vary the terms of this warranty.

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